

Extracting Vitals from Smartphone Accelerometer

I set out to detect heart rate and respiration using only a smartphone accelerometer. I conducted controlled breathing and spontaneous breathing experiments (2s, 4s, 8s cycles).




Data from an accelerometer is messy

1

Raw Z-axis Data


The raw signal is dominated by gravity, sitting at 9.8.

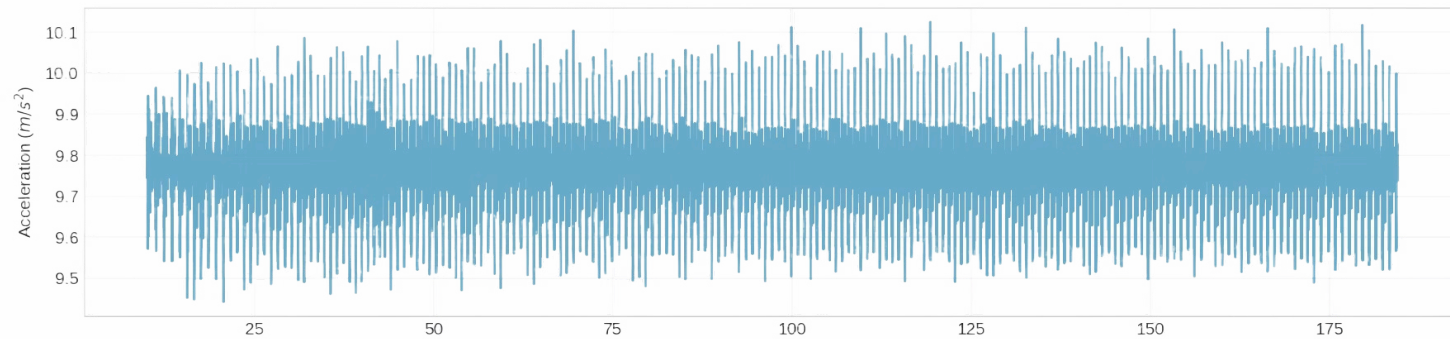


2

Intermediate Output: Raw PCA Signal (Pre-Filtering)


1. Detrend (Remove Gravity).
2. PCA (Combine Axes into Max Variance direction).
3. Bandpass Filter (Isolate Frequencies).

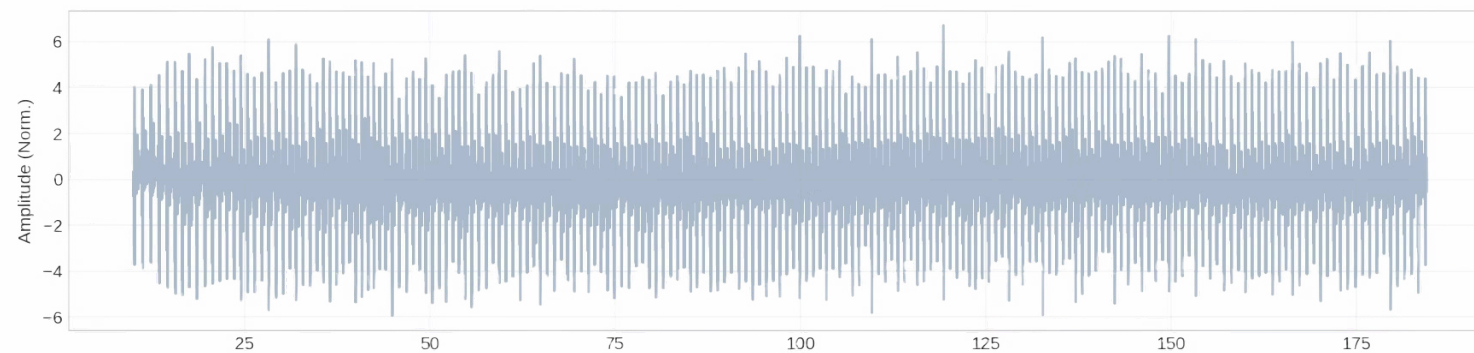

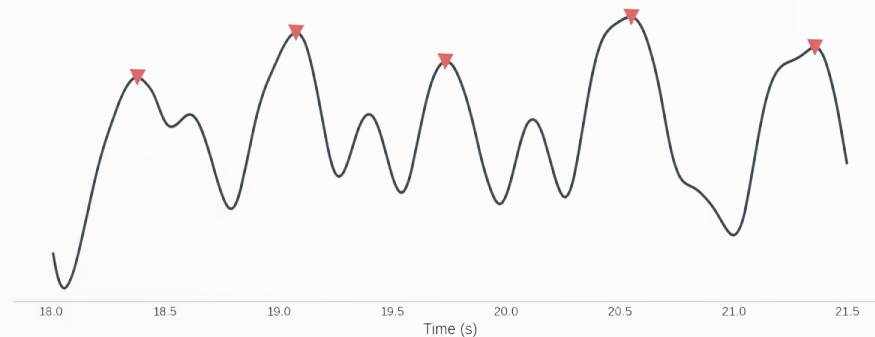


The Standard Pipeline (Time Domain)

- **The Goal:** Identify local maxima (peaks) as heartbeats.
- **The Math:** Calculate the average time difference between all valid peaks to derive a single Average Heart Rate for the recording.
- **The Assumption:** "If we filter the signal well, every peak is a heartbeat."
- **The Reality:** "Unfortunately, real-world data isn't this clean..."

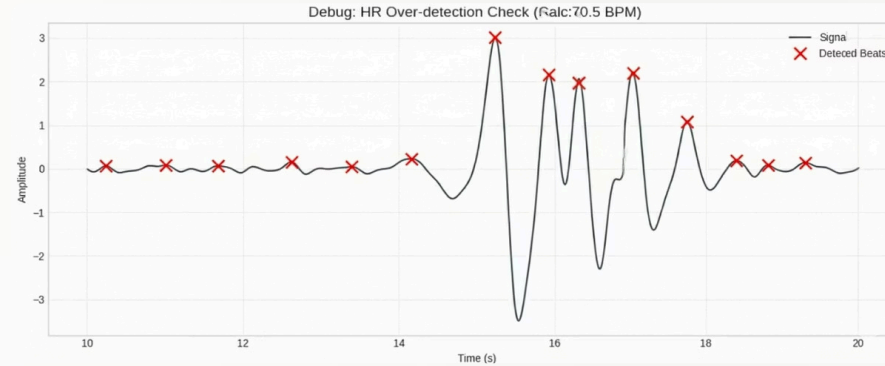


The Limitations of Static Thresholds

- **Method:** Standard peak detection with fixed amplitude.
- **Failure:** Motion artifacts are miscounted as heartbeats.
- **Result:** Over-estimation of Heart Rate.

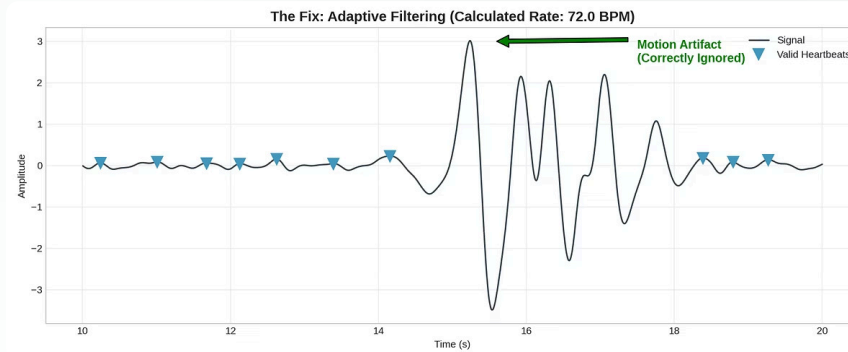


The Fix: Dynamic Thresholding

- **Physiological Constraints:** I implemented a **Minimum Peak Distance**. The algorithm now "knows" that a heart cannot physically beat twice in 0.2 seconds, preventing double-counting.
- **Median-Based Rejection:** I switched to **Median** filtering to calculate the "typical" peak height, ignoring massive motion spikes.

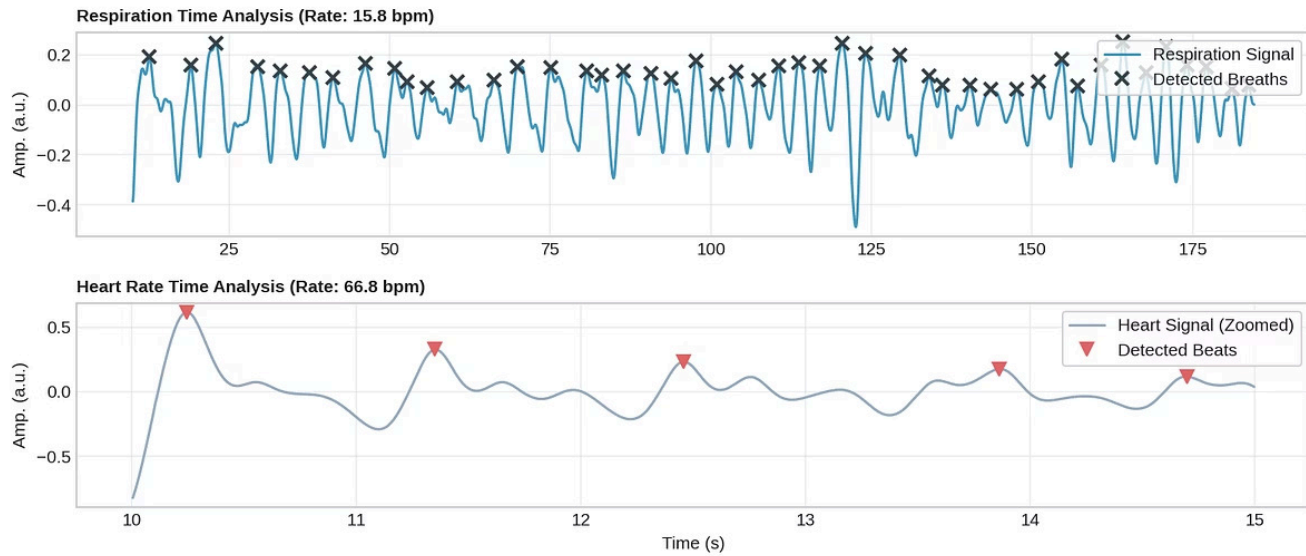


Final Vitals Extraction

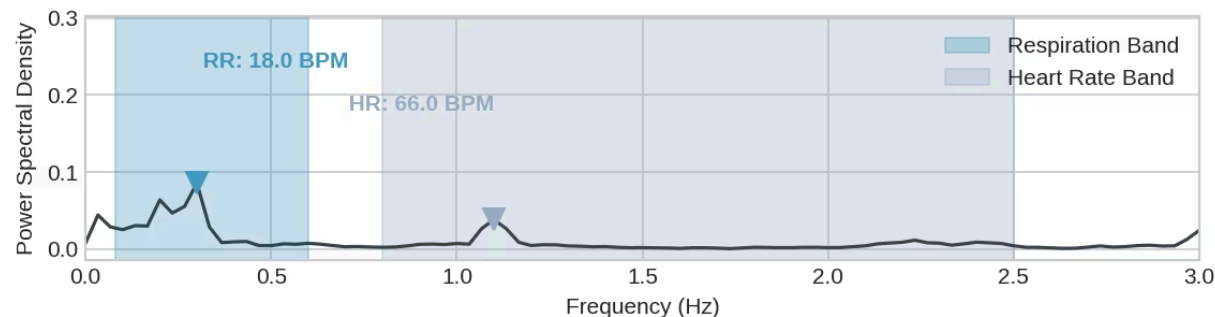

The Robust Solution: Welch's Periodogram

Switching to Frequency Domain analysis, specifically **Welch's** method, provided a breakthrough. This technique focuses on identifying dominant, repeating patterns across the entire signal rather than isolated spikes.

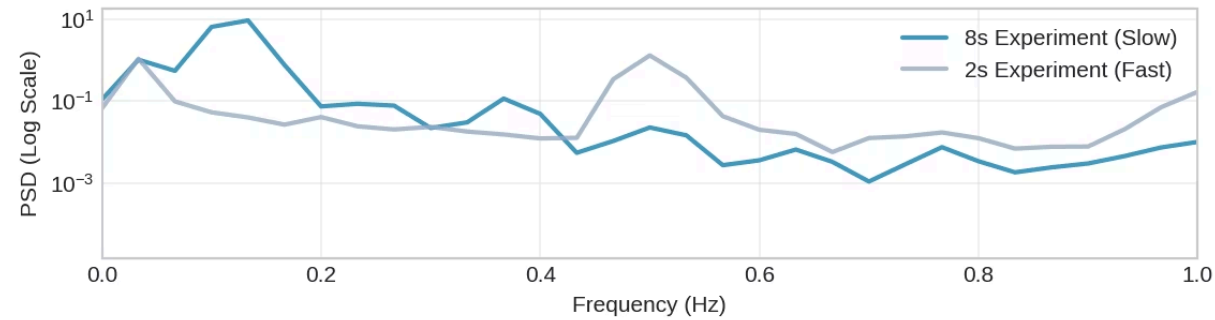

Proof of Robustness: Controlled Breathing Experiments

To prove this wasn't a fluke, I compared my 'Fast Breathing' (2s) experiment against my 'Slow Breathing' (8s) experiment.


Slow Breathing (8s cycle)

The **blue line** clearly shows a peak at approximately **0.125 Hz**, corresponding to the 8-second breathing cycle ($1/8$ Hz).

Fast Breathing (2s cycle)


The **grey line** demonstrates a distinct peak at **0.5 Hz**, aligning perfectly with the 2-second breathing cycle ($1/2$ Hz).

The **precise alignment** of these peaks with the expected frequencies proves that the algorithm accurately and reliably fulfils its purpose.

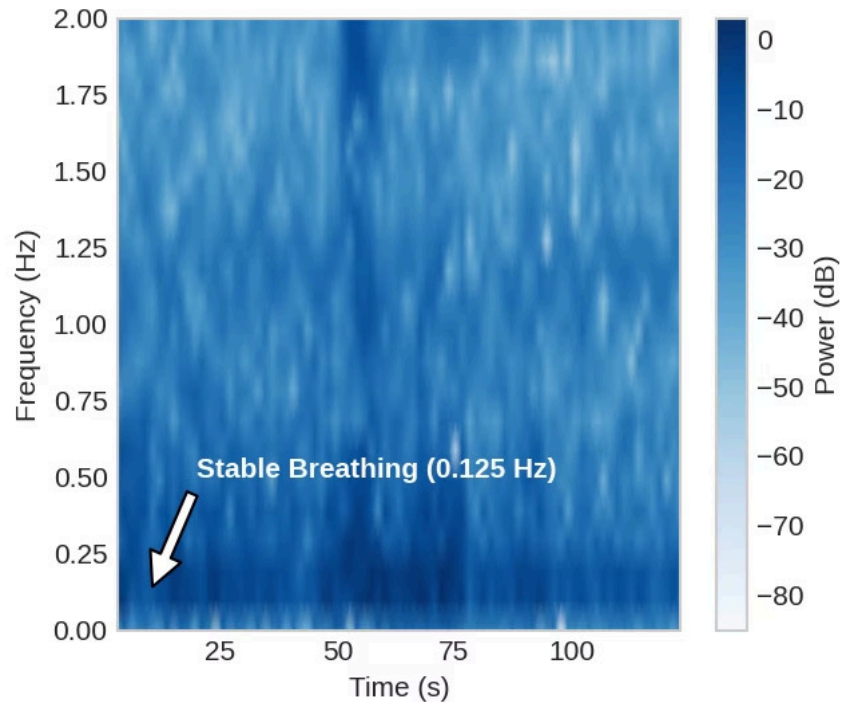

Visualising Stability: The Spectrogram

- **Unbroken Horizontal Line:** The spectrogram displays a consistent, continuous high-energy band across the entire 2-minute recording.
- **Experimental Control:** This visually confirms that the breathing rate **did not drift or fluctuate**, validating the quality of the controlled experiment.

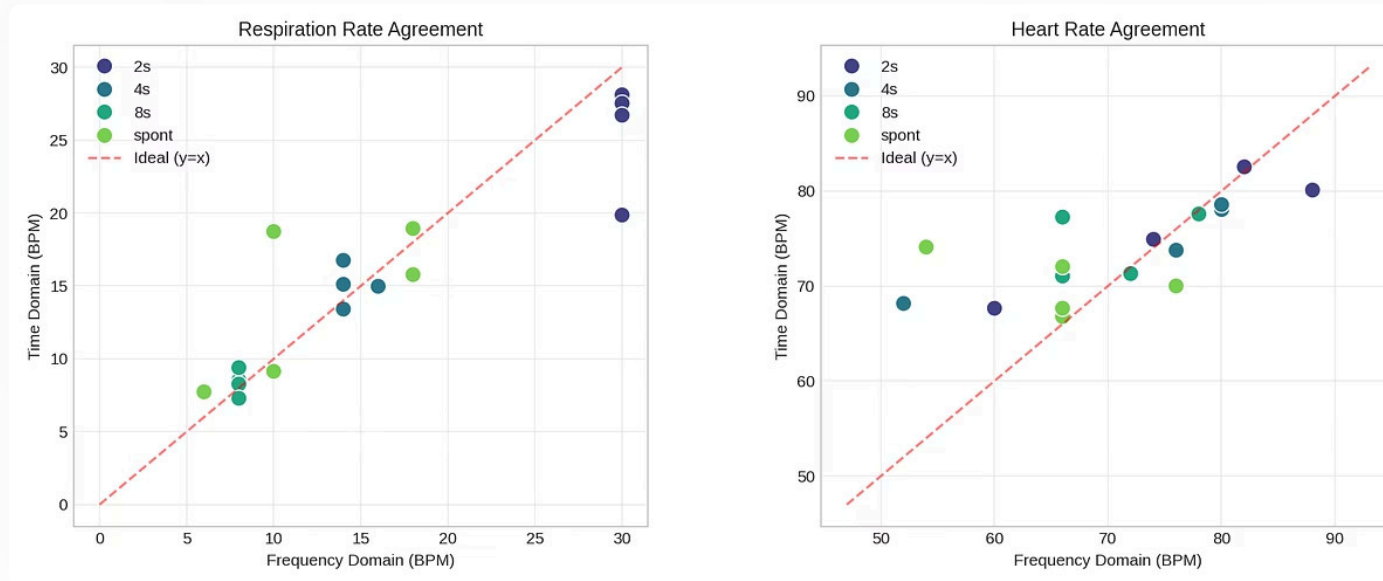

Validation: Method Agreement

- **General Trend:** The two methods show strong agreement, with data points clustering near the ideal 1:1 line for most experiments.
- **Heart Rate Bias:** Time Domain is slightly higher. Likely caused by residual false positives from noise artifacts increasing the count.
- **Respiration Bias:** Time Domain is lower during fast breathing. Shallow breaths sometimes fall below the detection threshold (false negatives). Even with dynamic thresholds some human intervention seems to be still necessary to adjust their parameters to ensure accuracy.

A decorative graphic on the left side of the slide, consisting of several overlapping, wavy, horizontal bands in various shades of blue, ranging from dark navy to light sky blue. The waves flow from the left towards the right, creating a sense of movement and depth.

Conclusion: Key Takeaways

- Vital signs *can* be extracted from noisy smartphone accelerometers, but it requires moving beyond naive assumptions.
- **Frequency Domain** is excellent for average rates (Stable)
- **Time Domain** is viable *if* adaptive thresholds are used.